




MANUAL WHOLE SLIDE IMAGING

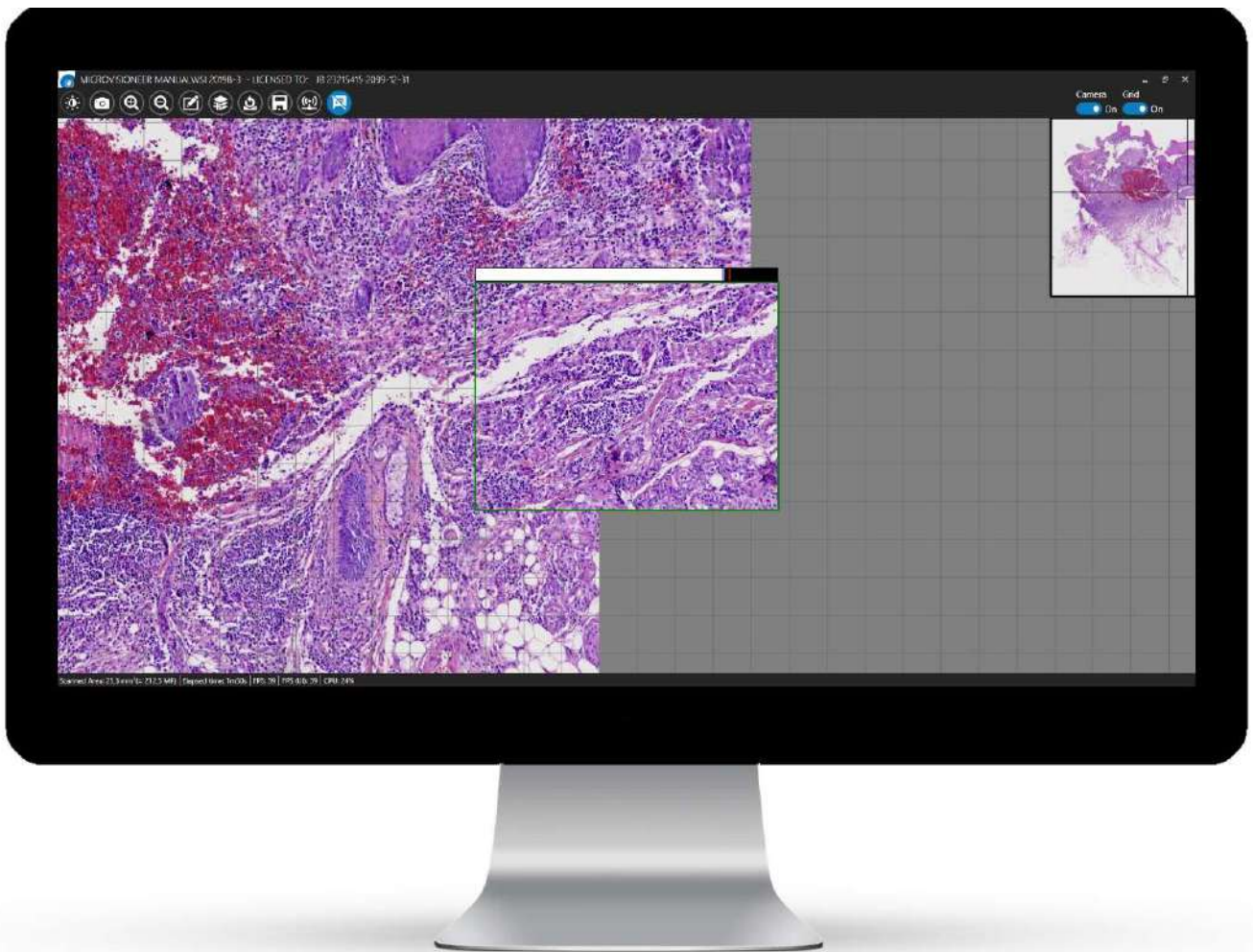

Microvisioneer mvSlide
The easy slide scanning solution

<https://www.microvisioneer.com/mvslide>



The Microvisioneer mvSlide Scan Software

Upgrade your microscope to a manual slide scanner

Microvisioneer mvSlide is an affordable and high-quality manual scanning software allowing to upgrade existing microscopes to manual whole slide scanners. The exceptionally versatile and innovative solution enables customers in 60+ countries worldwide to create virtual whole slide images (WSI) for various applications. The mvSlide system is affordable, easy to use and readily available, and allows smooth image stitching at high quality, similar to that of automated slide scanners. Almost any type of slide can be digitized, no matter how exceptional.



Microvisioneer offer an innovative software and personalised support package allowing users to convert digital images directly into whole slide images (WSI) of a quality comparable to commercial whole slide scanners.

Jeremy Parry,
Managing Director
DLS Pathology Services



The manual scanning process


- 1** Upgrade your microscope to a manual slide scanner with a suitable camera and the Microvisioneer mvSlide scan software. Camera mounting and software installation is extremely easy. A simple microscope upgrade – the ticket to virtual microscopy.
- 2** Create whole slide images by manually moving the stage. It has never been easier. Any desired objective can be selected. Scanning is intuitive and smooth. Real-time stitching of the panoramic image is seamless.


- 3** View, analyze, archive or share virtual slides. The virtual slides are compatible with all major available viewers and image analysis tools. For sharing the images, simply use cloud-based systems or the Microvisioneer InstaSlide software.



The Digital Workflow with mvSlide

The easy transition to fully digital workflows

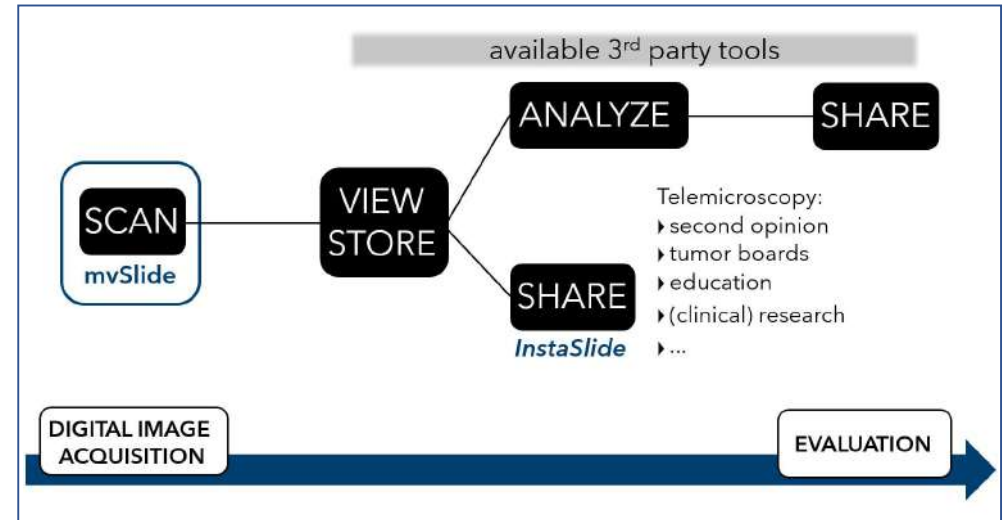

The Microvisioneer manual scanning solution can be readily integrated into any environment. The file formats of the whole slide images are standard and open industry file formats which are compatible with all major existing digital microscopy tools, such as viewers, image analysis software or cloud platforms. Thus, no major capital investment is necessary and a full and cost-effective telemicroscopy workflow can be established.

Microvisioneer mvSlide is an extremely powerful technique as it is easy to access and therefore significantly facilitates a transition to digital workflows.

For quick and straight-forward slide sharing, Microvisioneer also offer a telemicroscopy solution, *InstaSlide*, that can be integrated into mvSlide for fast communication with remote colleagues.

Benefits and Opportunities

Good reasons for manual scanning with Microvisioneer mvSlide




▶ Time-saving:

Hardware items and software are immediately available.
Manual scanning of small specimens is often faster than automated scanning.
Compensation of potential outage of automated scanners.
Scan right away on your desk!

▶ Affordable:

Upgrading of microscope with the mvSlide software is inexpensive.
It is a long-lasting solution without expensive maintenance.

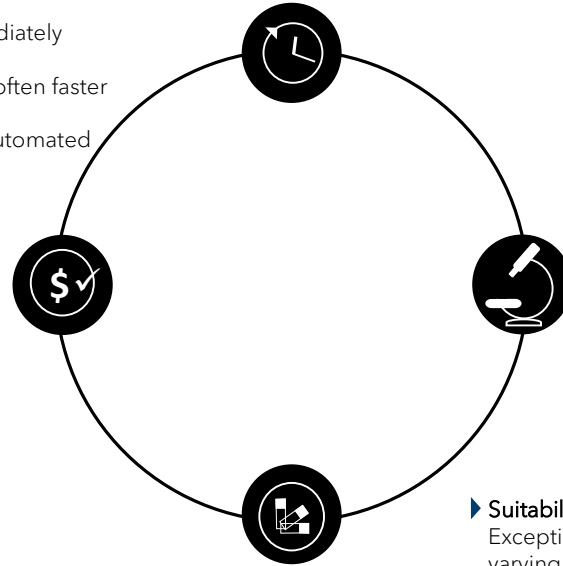

▶ Maintains traditional microscope features:

Any desired objective can be selected.
Slides can still be assessed through the eyepiece.

▶ Suitability for exceptional slides:

Exceptionally thick slides or slides with varying focus plane can be easily scanned due to the flexibility of the microscope.

mvSlide opens up all the opportunities of digital microscopy



Image Analysis

Perform representative and reproducible analyses.
Automate image analysis.



Documentation & Publication

Save your digital slides, annotations and comments.
Preserve the original staining.
Create publication-ready high-quality images.



Education

Scan your own teaching slides and move away from expensive glass slides.
All students see exactly the same virtual slide, remotely and simultaneously.




Collaboration


Share digital slides with colleagues worldwide.
Work remotely and save time and costs.

Applications


You decide what you scan

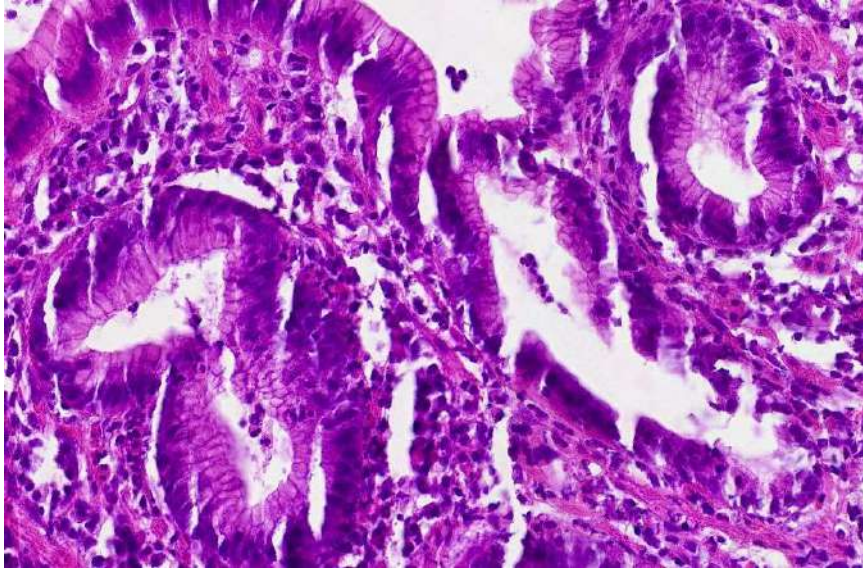


brightfield
human
intestine,
H&E stain

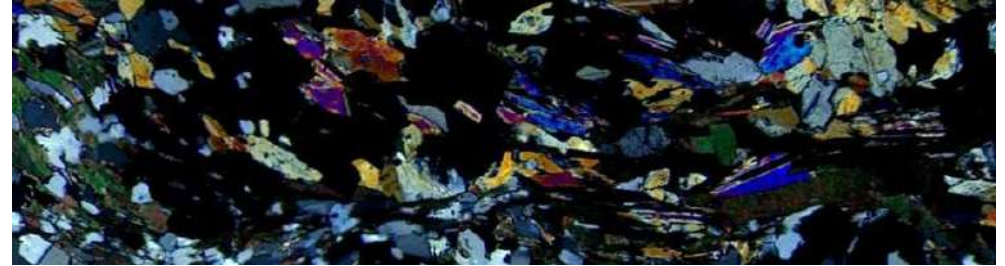

polarized light
rock thin section

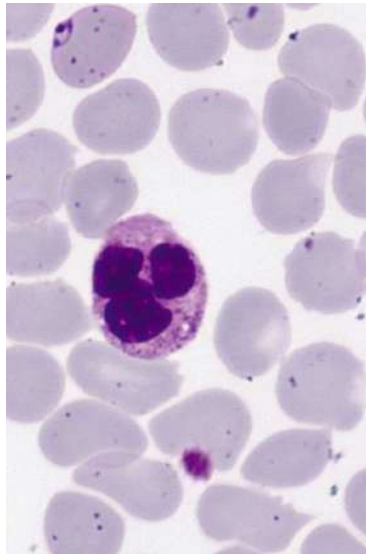


**oil
immersion**
blood
smear,
Wright's
stain

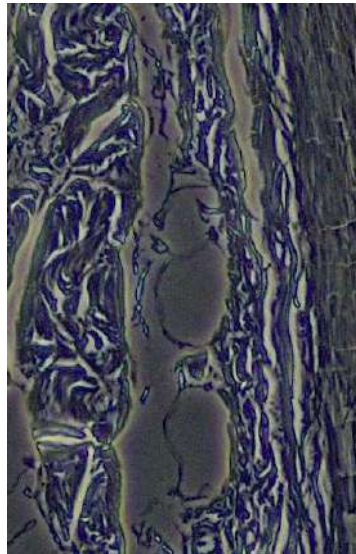

**phase
contrast**
skin of
pig
embryo

fluorescence
FluoCells™
prepared
slide #3
(mouse
kidney
section with
Alexa Fluor™
488 WGA,
AlexaFluor™
568
Phalloidin,
and DAPI)

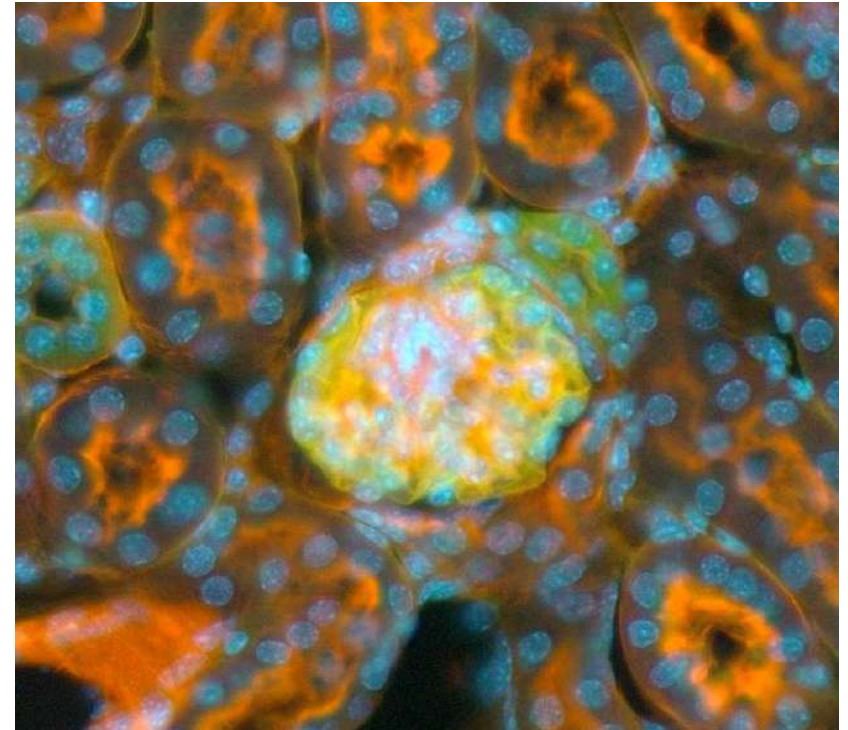

Specifications

Key technical features of Microvisioneer mvSlide

Supported microscopes	Trinocular with PLAN lenses or better
Supported scan modes	Brightfield, darkfield, phase contrast, polarization, oil immersion, and fluorescence
Supported cameras	Cameras with global shutter sensors or CCD sensors in progressive scan mode
Output file formats	SVS, TIFF, JPG, BMP, PNG
Operating system	Win7 or higher, 64bit
Shading compensation	Active compensation of shading / vignetting
Flicker compensation	Active compensation for older microscopes with flickering light sources
Freeze feature	Combination of multiple scans of different tissues / stains into one virtual slide
Snapshot feature	Shortcut to automatically save still images with or without scale bar while scanning. Can also be triggered via USB foot switch.
Presets	Store and load settings for different objective lens magnifications and stains
Realignment	Semi-automatic, automatic
Barcode reading for high-throughput scanning	Integration of barcode reading for automated file name creation including regular expressions and time stamps
Label and Macro Add-On	Add-on USB camera to capture whole slide and label, inclusion of the image in the scan
Pen annotation feature	Direct integration of basic annotations into the scan
White balance	Automatic and manual
Exposure adjustment	Automatic and manual
Possibility of integrating Microvisioneer InstaSlide and mvLive	Microvisioneer's solutions for sharing scans and for broadcasting the live camera image can be directly integrated

Microvisioneer mvSlide is for research use only.


Microvisioneer Product Range

Find your solution



Microvisioneer mvSlide Brightfield Edition

Manual scanning for brightfield applications, oil immersion, polarized light and phase contrast. Suitable for histology, immunohistochemistry, cytology, material science and many more.

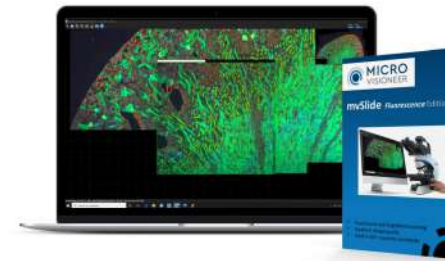

Microvisioneer mvSlide Fluorescence Edition

Manual scanning for *brightfield and fluorescence* applications. Both single channel and multi channel fluorescence imaging is possible.



Microvisioneer InstaSlide Slide Sharing Software

Immediate sharing of whole slide images via the internet with anyone, anywhere. A Microvisioneer telemicroscopy solution.



Microvisioneer mvLive Broadcasting Software

Live sharing of the camera image with remote colleagues. Best image quality without compression artefacts.

Microvisioneer GmbH

Germany

VAT ID: DE366777766



Contact us to discuss your opportunities or to book an online demo.

Email info@microvisioneer.com

Phone +49(0) 151 41605950

<https://www.microvisioneer.com>

